

ISOCAL - 6 Range Hyperion • Drago

- 65 x 160mm Diameter Calibration Volume
- Portable Liquid Bath
- Convertible for Dry Block Operation
- Calibrate all Sensor types

The Hyperion and Drago have a large calibration volume, 65mm x 160mm deep, which makes them ideal to use as portable liquid baths. Stirred liquid baths are suitable for temperature sensors of all types, sizes and shapes.

Liquid Baths can provide smaller calibration uncertainties than dry blocks and, when used with suitable reference thermometers, accuracies of up to 0.005°C can be achieved. The sensors can be placed directly into the stirred liquid thus avoiding the need for specially drilled blocks.

An ideal combination would be to use one of the Isotech TTI temperature indicators and semi standard probe, for example the 935-14-16 Semi Standard.

Alternatively, the SITE model includes an inbuilt temperature indicator and can be provided with a reference probe and calibration certificate thus giving a complete traceable solution in a single portable unit.

These models are part of the Isocal-6 family and with a reference probe can be used with different accessories for Dry Block, Infrared, Surface Calibration and even with ITS-90 Fixed Point Cells for uncertainties to 0.001°C.

In the optional Dry Block Mode, the large 65mm diameter block allows for the calibration of either larger probes or for calibrating many sensors simultaneously over the range -25°C to +250°C. Fully compatible with the leading I-Cal Easy automatic calibration software you can also benefit from the combination of award winning flexibility and large calibration capacity.



<http://www.isotech.co.uk/industrial/>



Specification

Model	936 Hyperion ^{PLUS}	934 Drago ^{PLUS}
Temperature Range	-25°C to 140°C (1)	30°C to 250°C (2)
Stability	Dry Block: ±0.03°C, Blackbody: ±0.3°C, Surface Sensor: ±0.5°C, Liquid Bath: ±0.025°C, ITS-90 Cells: ±0.0005°C, Ice Bath: ±0.001°C (not Calisto)	
Accuracy (3)	0.15°C	0.15°C
Uniformity - Between Wells Dry Block Mode (Radial)	<0.008°C	<0.008°C
Uniformity - Radial Liquid Bath Mode	<0.009°C	<0.007°C at 250°C
Uniformity - Lower 40mm (Axial) Dry Block Mode	<0.040°C	<0.040°C
Uniformity - Lower 40mm (Axial) As Liquid Bath	<0.011°C	<0.013°C at 250°C
Heating Time	-20°C to 140°C: 40 Mins	30°C to 250°C: 40 Mins
Cooling Time	140°C to 20°C: 90 Mins 20°C to -25°C: 80 Mins	250°C to 30°C: 90 Mins
Insert Diameter	65mm	
Immersion Depth	160mm	
Insert Types	Standard 8 x 8mm + 2 x 4.5mm, Undrilled or Custom	
PC Interface	Included - Supplied with PC Cable and Software	
Power	200 Watts	1000 Watts
Voltage	115Vac or 230 Vac 50/60Hz	
Dimensions	302 x 176 x 262 mm	
Weight	12kg	8kg

(1) In ambient of 20°C: Minimum Temperature is 45°C Below Ambient, Absolute Minimum -35°C

(2) In ambient of 20°C

(3) Dry Block Mode only: Comparing 4.5mm Well to Controller Display Value.

A three point traceable calibration certificate is included. UKAS calibration options available, contact Isotech for advice and options.

Features (Basic & Site)	936 Hyperion ^{PLUS}	934 Drago ^{PLUS}
Dry Block	✓	✓
Stirred Liquid Bath Option	✓	✓
Stirred Ice Bath Operation	✓	✓
Surface Sensor Option	✓	✓
Infrared Calibration Option	✓	✓
ITS-90 Fixed Point Cells	Water, Gallium	Gallium
Additional 8mm Pre-heat Pocket	✓	✓
Configurable Units: °C, °F and K	✓	✓
Supply Voltage Power Correction with Digital Filtering	Provides High Stability protecting against noise and supply voltage variation	
Additional Features (Site)		
Independant Temperature indicator	✓	✓
Universal Input Types PT100	✓	✓
Thermocouples Types K,N,R,S,L,PL2,T,J,E	✓	✓
Linear Process Inputs Including 4-20 mA	✓	✓
Stand Alone Thermostat Testing	✓	✓
Thermostat Testing With PC	✓	✓
Five Point Digital Probe Matching	✓	✓
Configurable Units: °C, °F and K	✓	✓

Isocal-6 Ultimate Flexibility - Hyperion & Drago

Calibrate all sensor types - Thermocouples, PRT's, Thermistors, Thermostats, Infrared, Surface Sensors...



1 Dry Block Mode with Inserts
936-06-01a Standard Insert is:
 8 x 8mm + 2 x 4.5mm all 157mm Deep.
 All Inserts have a 4mm tapped hole to
 suit supplied extractor tool.

Alternative Inserts

- 936-06-01b** Blank Insert
- 936-06-01c** Special Insert.
*Contact Isotech with
 your requirements.*




**5 Surface Sensor Calibration with
 Surface Sensor Kit**
936-06-04 Includes an Insert and
 an angled platinum
 resistance thermometer.




**2 Stirred Liquid Mode
 with Liquid Container Kit**
 Allows liquid bath use, includes
 container, magnetic stirrer, probe guide
 and sealing cap. **936-06-02**

**3 Stirred Ice Bath Mode with Liquid
 Container Kit**
 Uses same liquid kit to provide 0°C
 reference as a stirred ice bath.

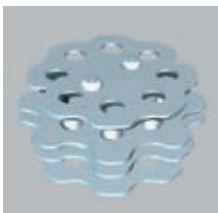


**Additional Accessories for use with
 the above kit**
Thermometer Support Kit
936-06-08 Supports up to eight
 thermometers into liquid.
 Suits probes 5mm - 8mm
 in diameter.

- Oils**
- 936-06-07** C10 Oil -35°C – 140°C 1L
 - 580-06-09** C20 Oil 20°C – 200°C 1L
 - 915/09** VH Oil 150°C – 250°C 1L

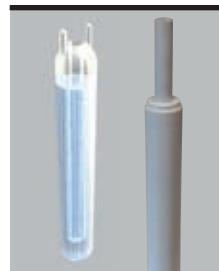


6 ITS-90 Fixed Point Cells
B8 Water Triple Point Cell
 (Hyperion)
17401 Slim Gallium Slim Cell
936-06-09 Cell Holder Assembly

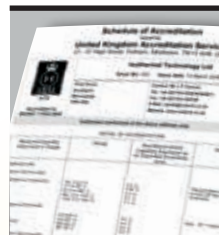

UKAS Calibration
 UKAS Calibration available to order,
 legally traceable in more than 70
 countries.



Standard Probe
935-14-61/DB Platinum Resistance
 Thermometer.
 4mm diameter.

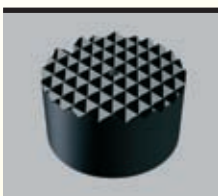

**4 Infrared Calibration Mode with
 Blackbody Target**
936-06-03 Use optional Probe
 936-14-61DB.




Carrying Case
931-22-65 Sturdy case
 accommodates the unit
 with room for
 accessories


How To Order
 Specify Model, Basic or Site, Supply Voltage, Accessories and
 if UKAS Calibration is required.

Isocal-6 Performance and Use

Hyperion^{PLUS} Performance - Dry Block

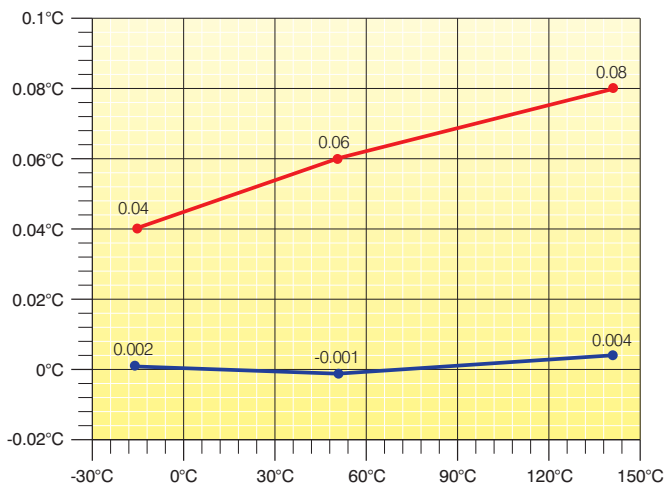

Drago^{PLUS} Performance - Dry Block

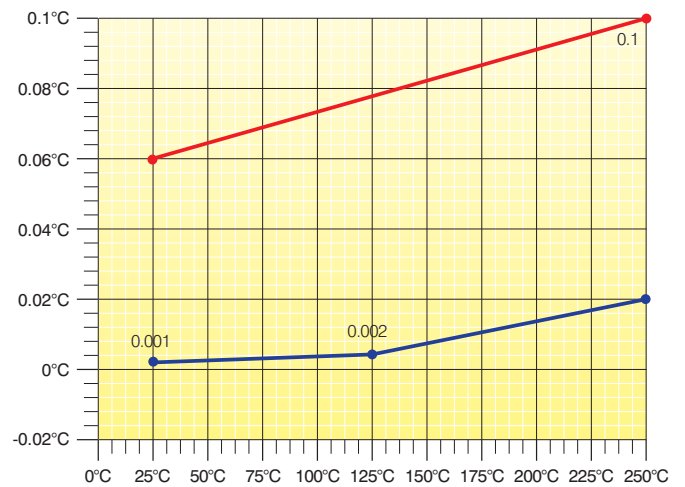

- Audit Calibration (Similar Sensors) S model with UKAS option
- Radial Homogeneity

See Evaluation Reports for full details
<http://www.isotech.co.uk/refer.html>


Alternative Methods of Calibrating with an Isocal-6



A Basic Dry Block Calibrator

For quick and easy Testing and Calibration. Relies on using the temperature controller value, "controller accuracy".


Using a standard thermometer (not shown) compensates for temperature gradients and loading errors giving best performance.

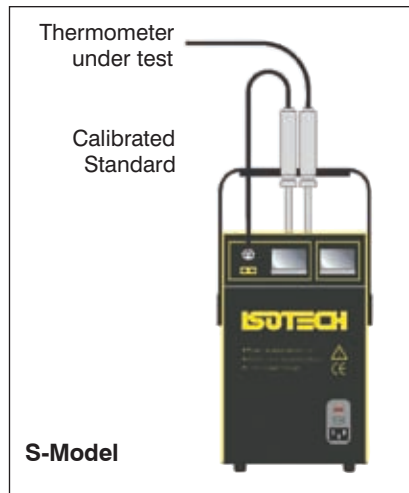

An ISO 9000 Calibration System

Thermometer under test is compared to a calibrated standard and inbuilt independent temperature indicator.

For best practice calibration with established traceability and uncertainty.

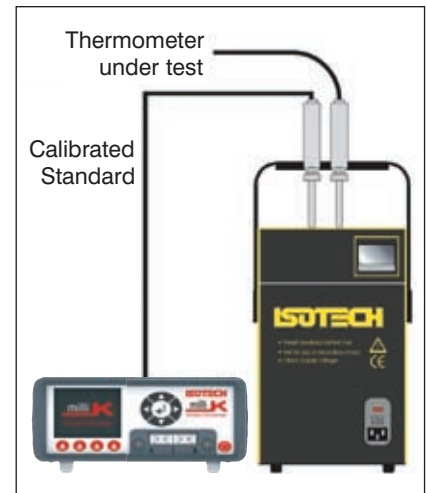

Using an External Indicator

Thermometer under test is compared to a calibrated standard and external temperature indicator.

For best practice calibration with established traceability and uncertainty.

One indicator may be used with several Isocal-6 or similar.

Calibration is needed only for the probe and indicator.

For liquid bath, surface sensor or black body always use a standard thermometer connected to the S-model or stand alone indicator as above.



UKAS Calibration available for these systems - *International Traceability - Best Practice*