

## Flammable and combustible liquid storage

- Upscale design and sized for safe storage of flammable liquids in laboratories
- Thermal resistant, 1<sup>1</sup>/2" (38mm) wood construction with attractive stainless silver color or Light Neutral powder paint finish
- Removable 12" (305mm) square rear panel allows access to plumbing or electrical connections
- Adjustable leveling feet are accessible from inside the cabinet for added convenience
- Self closing doors with built-in fusible links will release in the event of a fire. Ball lock mechanism holds doors open for easy access
- Door access panel and vent openings are lined with intumescent seal designed to swell and fill the gap in the event of a fire
- Each door individually locks with matching keys
- Steel components are bonded by a common strip all connected to a ground/bond terminal located on the back of the cabinet
- Supplied with one adjustable shelf 1.22" (33mm) deep
- Generous sump watertight to 4" (102mm) deep
- Haz-Alert<sup>™</sup> reflective warning label
- Two diagonally opposed 2" (51mm) NPT vents on the back, shipped plugged



## U.L. Listed Safety Cabinets For LABORATORIES

These double-door steel-clad cabinets feature 1<sup>1</sup>/2" (38mm) thick thermal-resistant wood construction and are built to comply with stringent Underwriters Laboratories (UL) and FM Global (FM) standards. Equipped with self-closing, independently lockable doors with fusible links and an intumescent lining around all openings, these cabinets offer the optimum protection under fire conditions. A generous sump captures hazardous spills. Vent capability on the back of cabinet and a rear access panel for utilities offer added convenience in laboratories or undercounter settings.

> For added protection, intumescent seal fills any gaps under fire conditions.

High visibility, reflective Haz-Alert<sup>™</sup> warning label.










Convenient, removable rear access panel for utilities. Two 2" (51mm) flame arresting vents for optional venting.

Steel shelf and sump are bonded by a common strip inside the cabinet and connected to a ground terminal on back of cabinet. This allows for safe dispensing inside the cabinet.

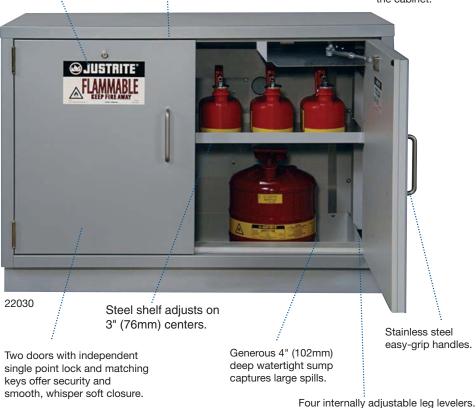

## U.L. Listed Laboratory Safety Cabinets New Items

	<b>Dimensions - inches</b>	Ship Wt.	Silver	Lt. Neutral
Capacity	HXWXD (mm)	Lbs./Kgs	Model No.	Model No.
13 Gallon/49L	35 <sup>3</sup> /4" x 30" x 21 <sup>5</sup> /8" (908 x 762 x 549)	247/112	22010	22011
17 Gallon/64L	35 <sup>3</sup> /4" x 36" x 21 <sup>5</sup> /8" (908 x 914 x 549)	293/133	22020	22021
25 Gallon/95L	35 <sup>3</sup> /4" x 48" x 21 <sup>5</sup> /8" (908 x 1219 x 549)	375/170	22030	22031



