

F-103 Gen-5 Bench Fume Hood

F-103 Isolator Gen-5 Bench Fume Hood - Framed Vertical Sash

ISOLATOR Gen-5 Bench Fume Hood - This general purpose fume hood is designed to meet most laboratory Constant Air Volume (CAV) or Variable Air Volume (VAV) requirements and supplied with and automatic compensating upper by-pass. For VAV, use option MOD-0004 for a restricted by-pass position of observation glass. Designed to mount on a 31" deep counter top. The ISOLATOR Gen-5 bench mounted fume hood is supplied with the following standard features:

Counter Balance Systems - Chain and sprocket system delivers the easiest and most reliable sash operation available with an exceptional product life-cycle ISOLATOR drive system provides economical and reliable operation; tested to over 1 million cycles with no failure.

Sash Design - 1/4" thick tempered safety glass is provided with an epoxy urethane coated steel frame. Optional laminated safety glass or polycarbonate available for special applications. Hood has a framed vertical sash with a 42" exterior sightline.

Lower Airfoil - Aerodynamic flush sill design allows air to enter the fume hood even when the sash is closed ensuring efficient fume exhaust. The flip up safety trough cover allows electrical wiring or tubing into the fume hood while still permitting full closure of the sash.


Stainless Steel Type 304 Exhaust Collar - Round collar which allows for direct connection to typical exhaust ducts to reduce duct transition costs, minimize static pressure losses and exhaust noise levels.

Electrical - Two UL/CSA approved duplex receptacles provided for 120 volt service, one on each corner post. UL/CSA approved LED light fixture and switch provided. Explosion proof fixtures are available options.

Plumbing - Both corner posts are pre-punched to accept a maximum of five plumbing fittings per post. Factory pre-plumbing is available as well as plumbing fixtures from a variety of manufacturers to meet most plumbing needs.

Access Panels - Interior overlay style access panels provide convenient access and prevent accidental removal. Exterior side panels are also removable for ease of access to plumbing and electrical service fixtures.

Standards - UL 1805 Classified and tested in accordance with ASHRAE 110. Test results available upon request.


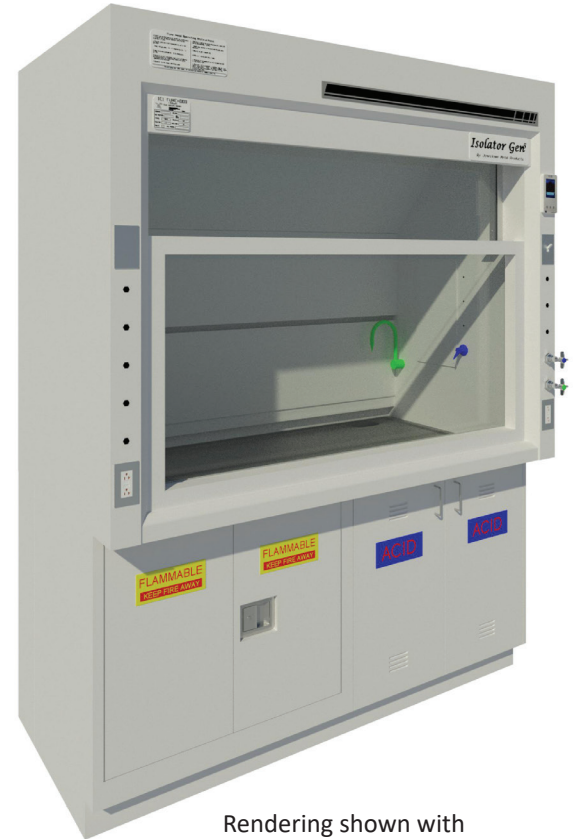


Rendering shown with optional features



Catalog Number	Width	Clear Interior Depth	Duct Diameter	28" Vertical Open @ 100 FPM	18" Vertical Open @ 100 FPM
F-103-36G5	36"	24"	1:10"	506 CFM @ .14" SP	325 CFM @ .06" SP
F-103-36G5D30	36"	30"	1:10"	506 CFM @ .14" SP	325 CFM @ .06" SP
F-103-48G5	48"	24"	1:10"	739 CFM @ .29" SP	475 CFM @ .12" SP
F-103-48G5D30	48"	30"	1:10"	739 CFM @ .29" SP	475 CFM @ .12" SP
F-103-60G5	60"	24"	1:12"	972 CFM @ .27" SP	625 CFM @ .11" SP
F-103-60G5D30	60"	30"	1:12"	972 CFM @ .27" SP	625 CFM @ .11" SP
F-103-72G5	72"	24"	1:12"	1206 CFM @ .42" SP	775 CFM @ .17" SP
F-103-72G5D30	72"	30"	1:12"	1206 CFM @ .42" SP	775 CFM @ .17" SP
F-103-96G5	96"	24"	2:10"	1672 CFM @ .37" SP	1075 CFM @ .15" SP
F-103-96G5D30	96"	30"	2:10"	1672 CFM @ .37" SP	1075 CFM @ .15" SP
F-103-120G5	120"	24"	2:12"	2139 CFM @ .33" SP	1375 CFM @ .14" SP
F-103-120G5D30	120"	30"	2:12"	2139 CFM @ .33" SP	1375 CFM @ .14" SP
F-103-144G5	144"	24"	2:12"	2606 CFM @ .49" SP	1675 CFM @ .20" SP
F-103-144G5D30	144"	30"	2:12"	2606 CFM @ .49" SP	1675 CFM @ .20" SP

F-103 Gen-5 Bench Fume Hood



F-103-36G5 - 36" Wide

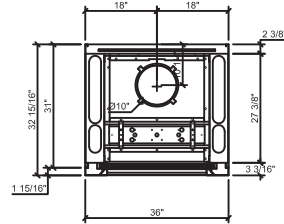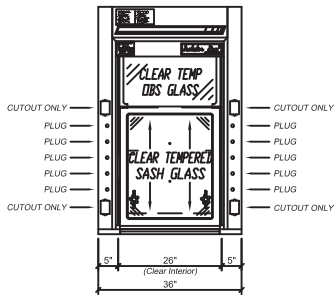

F-103-48G5 - 48" Wide

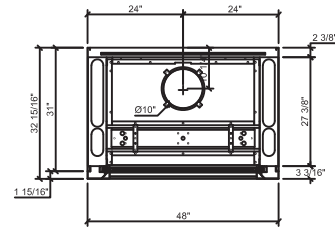
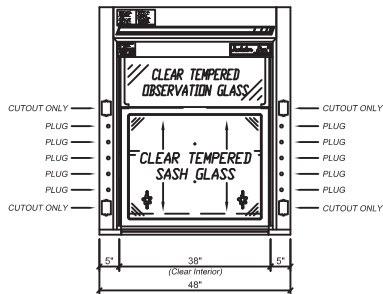


F-103-60G5 - 60" Wide

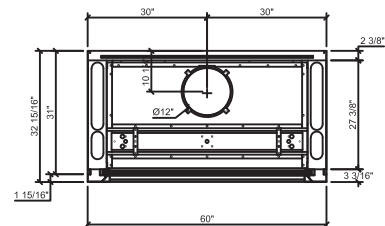
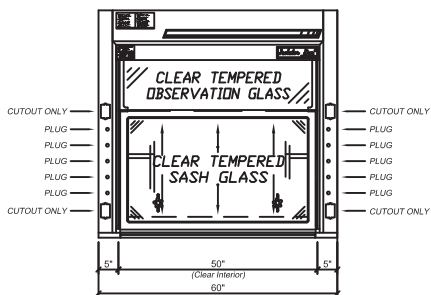


F-103 Gen-5 Bench Fume Hood

F-103-72G5 - 72" Wide

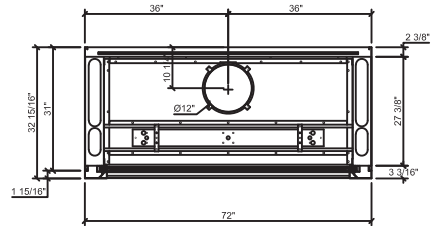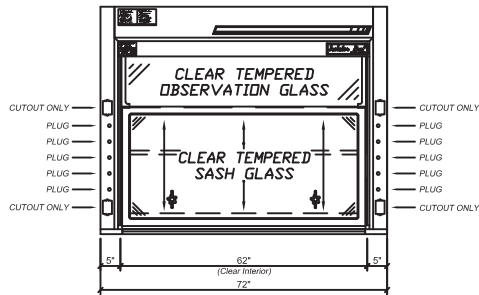

F-103-96G5 - 96" Wide

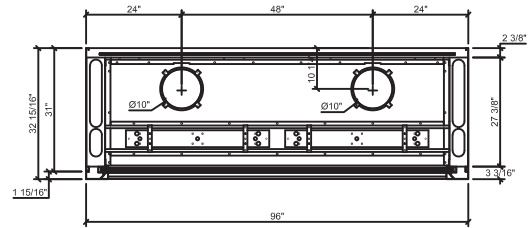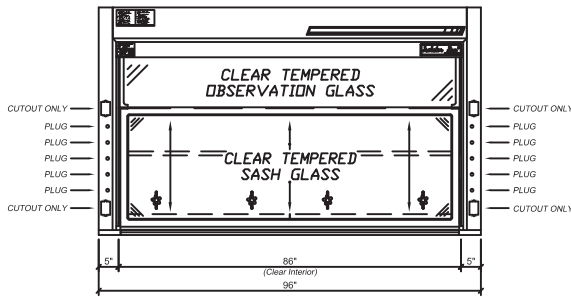

F-103-120G5 - 120" Wide

